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Hello QUA

```
from qm.qua import *
from configuration import qm

#####
# The QUA program #
#####
with program() as prog:
    with infinite_loop_():
        play("pulse", "rf1")

#####
# Execute the program #
#####
job = qm.execute(prog)
```

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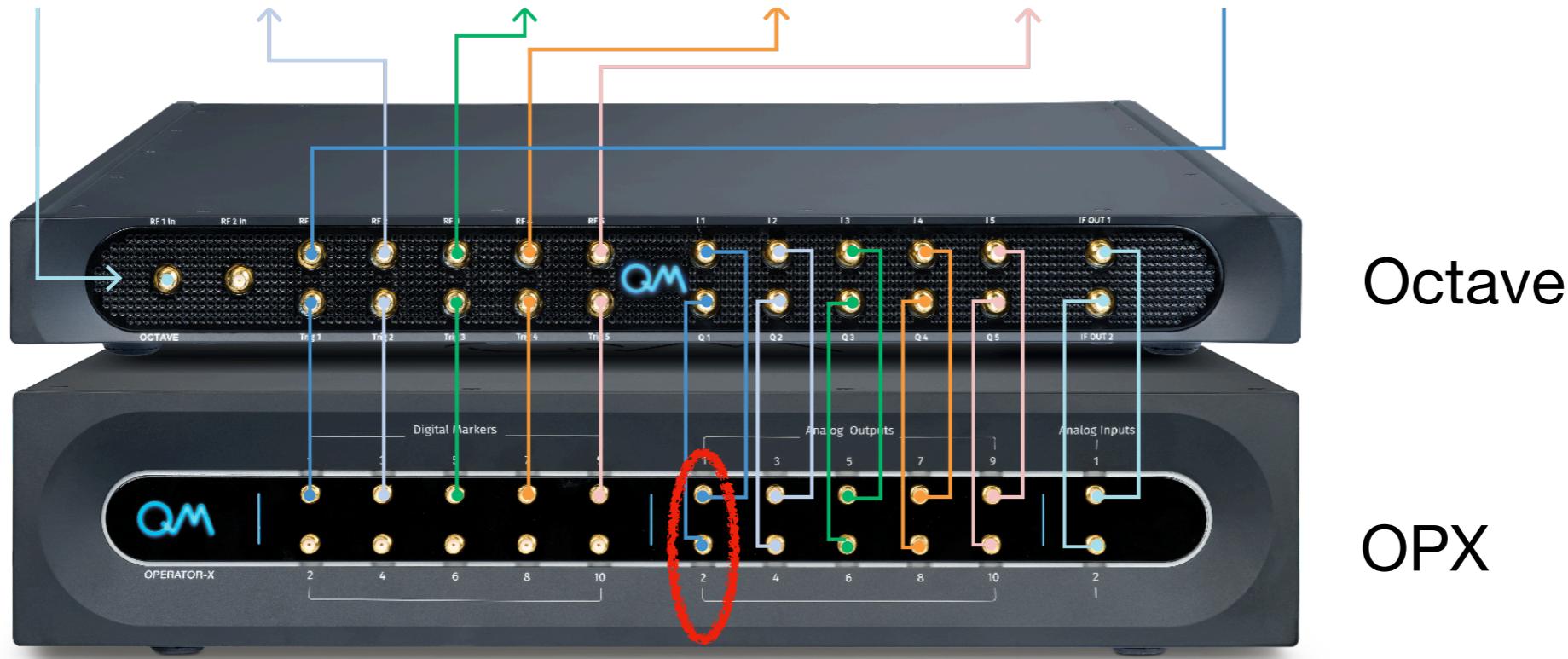
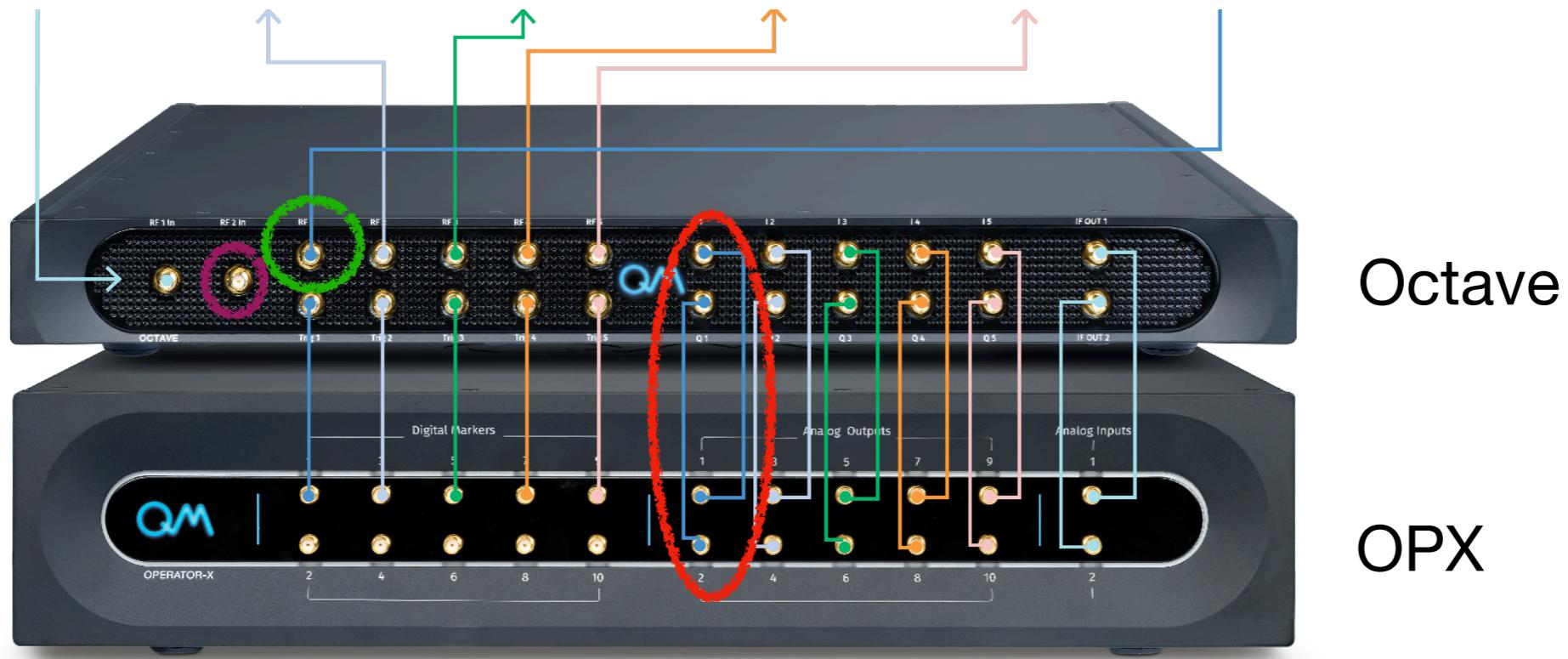
```
with program() as prog:
    I = declare(fixed) # QUA variable for the measured 'I' quadrature
    Q = declare(fixed) # QUA variable for the measured 'Q' quadrature
    I_st = declare_stream() # Stream for the 'I' quadrature
    Q_st = declare_stream() # Stream for the 'Q' quadrature

    with infinite_loop_():
        play("pulse", "rf1")

    with infinite_loop_():
        # Demodulate the signals to get the 'I' & 'Q' quadratures)
        measure(
            "readout",
            "spectro",
            None,
            dual_demod.full("cos", "sin", I),
            dual_demod.full("minus_sin", "cos", Q),
        )
        # Save the 'I' & 'Q' quadratures to their respective streams
        save(I, I_st)
        save(Q, Q_st)

    with stream_processing():
        # Cast the data into a 1D vector, and store the results on the OPX processor
        I_st.buffer(n_points).zip(Q_st.buffer(n_points)).save(`IQ`)
```

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
```
#####
# Live plot      #
#####
class myLivePlot(LivePlotWindow):
    def create_axes(self):
        # Create plot axes
        self.ax = self.canvas.figure.subplots()
        # Plot
        self.spectrum = self.ax.plot((spectro_L0 + frequencies)/1e6,np.ones(len(frequencies)))[0]
        self.ax.set_xlabel('Frequency (MHz)')
        self.ax.set_ylabel('Signal (dB)')
        self.ax.set_ylim(-120,-10)

    def polldata(self):
        # Fetch the raw ADC traces and convert them into Volts
        IQ = self.job.result_handles.get("IQ").fetch(1)
        if IQ is None:
            return
        I = IQ['value_0']
        Q = IQ['value_1']
        self.spectrum.set_ydata(10*np.log10(I**2+Q**2))
        self.canvas.draw()

#####
# Execute the program #
#####
job = qm.execute(prog)
qm.set_io1_value(1.0)
window = myLivePlot(job)
window.show()
```

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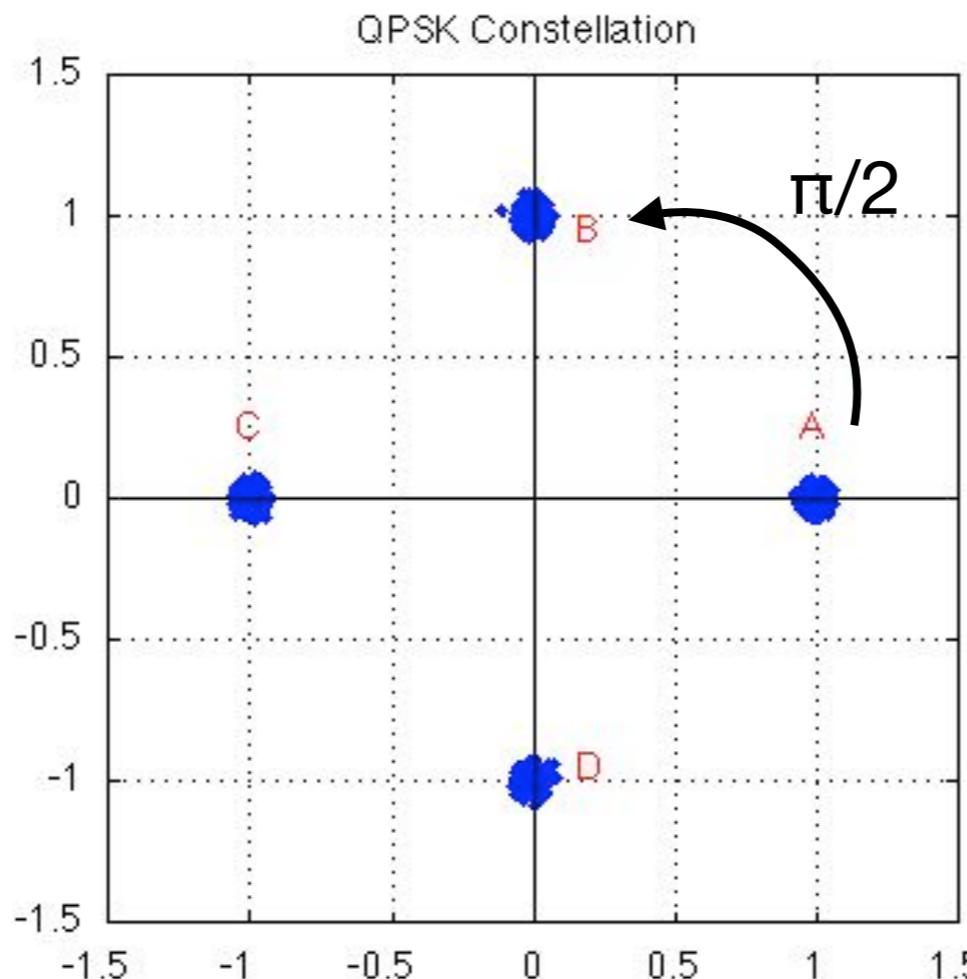```
msg=np.load("msg_1007.npz")
msgI=msg["x"]
msgQ=msg["y"]

with program() as prog:
    n = declare(int) # QUA variable for the averaging loop
    m = declare(int)
    Im = declare(fixed,value=msgI)
    Qm = declare(fixed,value=msgQ)
    I = declare(fixed) # QUA variable for the measured 'I' quadrature
    Q = declare(fixed) # QUA variable for the measured 'Q' quadrature
    I_st = declare_stream() # Stream for the 'I' quadrature
    Q_st = declare_stream() # Stream for the 'Q' quadrature

    with infinite_loop_():
        with for_(m, 0, m < len(msgI), m + 1):
            play("pulse" * amp(Im[m],0.,0.,Qm[m]), "rf1")
```

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